# PHILIPS

## RESPIRONICS

### DreamWear

#### Full Face Mask

# Instructions for use –

## 1 Safety Information

---

### **Intended Use**

This mask is intended to provide an interface for application of CPAP or bi-level therapy to patients. The mask is for single patient use in the home or multi-patient use in the hospital/institutional environment. The mask is to be used on patients (>66lbs/30kg) for whom CPAP therapy or bi-level therapy has been prescribed.

---

**⚠Caution:** U.S. Federal law restricts this device to sale by or on the order of a physician.

---

**🚫Note:** This mask is not made with natural rubber latex or DEHP.

---

### **Contraindications**

This mask is not for use on patients with the following conditions: recent eye surgery or dry eyes, hiatal hernia, excessive reflux, impaired cough reflex, and impaired cardiac sphincter function. This mask is not for use on patients who are dependent on mechanical ventilation for their life support. This mask is not for use on patients who are taking a prescription drug that induces vomiting, or on patients who are uncooperative, unresponsive or unable to remove the mask by themselves.

---



### **⚠Warnings**

- Use of this mask while the system is not on and operating may cause rebreathing of exhaled air. Rebreathing of exhaled air for longer than several minutes can in some circumstances lead to suffocation.

- Discontinue use and contact your healthcare professional if any of the following symptoms occur: skin redness, irritation, discomfort, blurred vision, or drying of the eyes.
- Discontinue use and consult your healthcare professional if any of the following symptoms occur: unusual chest discomfort, shortness of breath, severe headache, eye pain or eye infections.
- Discontinue use and consult your healthcare professional if any of the following symptoms occur: tooth, gum, or jaw soreness. Use of a mask may aggravate an existing dental condition.
- A minimum of 3 cm H<sub>2</sub>O (hPa) must be maintained when using this mask.
- Do not overtighten the headgear straps. Watch for signs of overtightening, such as excessive redness, sores or bulging skin around the edges of the mask. Loosen the headgear straps to alleviate symptoms.
- The headgear clips and mask cushion contain magnets. Contact your healthcare provider before you use this mask. Some medical devices may be affected by magnetic fields. The magnetic clips in this mask should be kept at least 2 in. (50 mm) away from any active medical device with special attention to implanted devices such as pacemakers, defibrillators and cochlear implants.
- Do not use in or near magnetic resonance imaging (MRI) equipment.
- Do not block or seal the holes and exhalation ports.


## 2 What's in the Package

---

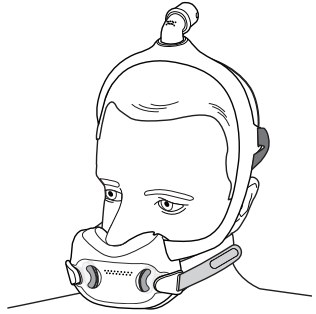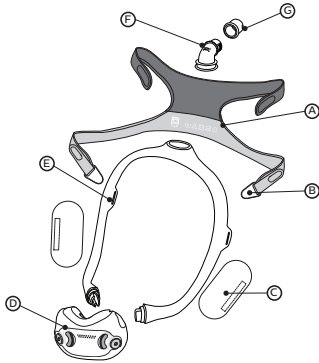
- Ⓐ Headgear
- Ⓑ Magnetic headgear clip
- Ⓒ Mask frame fabric sleeves (optional)
- Ⓓ Mask cushion (available in S, M, L, MW)
- Ⓔ Mask frame (available in small (SM), medium (MED), large (LG))
- Ⓕ Elbow with built-in exhalation (do not block)
- Ⓖ Tubing quick release (may remain in the CPAP hose)

Back of the Cushion

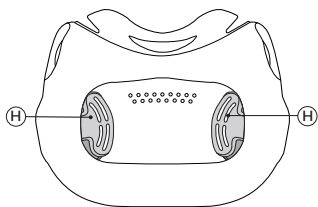Front of the Cushion

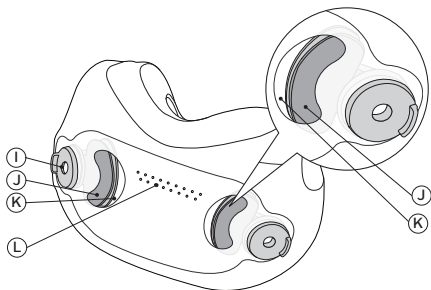- Ⓜ Guards
- Ⓜ Magnetic Clips
- Ⓜ Flaps
- Ⓜ Large Holes (Two) in the Cushion (do not block)
- Ⓜ Small Holes in the Cushion (do not block)


### 3 Before Use

#### DAILY BEFORE USE VERIFY FOR SAFETY

1. Look at the front of your mask cushion.
2. Find the two large holes.
3. Find the flaps. The flaps are clear, thin pieces that sit on the guards behind the holes.
4. Do not use if the two large holes are blocked shut by the flaps or anything else (for example, mucus). ①

**⚠Warning:** Do not block any of the mask holes.

**🚫Note:** The flaps should cover the two large holes when the device is on.
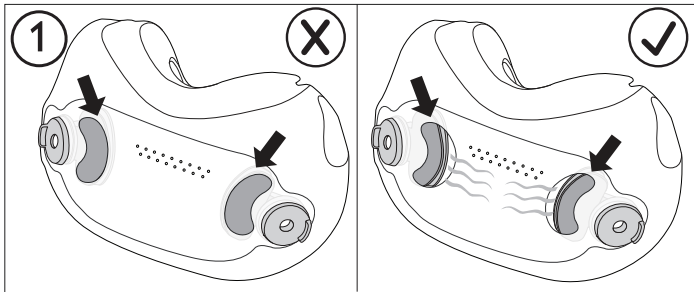
#### Before Use Read and Understand the Instructions Completely.

- Hand wash the entire mask before use (see **Mask Care** section).
- Wash your face. Do not use moisturizer or lotion on your hands or face.
- Inspect the entire mask daily. Discard and replace any damaged or worn parts.

## 4 Achieving the Right Fit

---

### Cushion Sizing

Four cushion sizes (S, M, L, MW) are available. Place the cushion sizing gauge under your nose (as shown).

**\* Tip:** For the best performance, use the smallest cushion that fits your nose.



1. Hold the gauge horizontally under your nose and select the cushion size based on where the outer edge of your nostrils and tip of your nose contact the gauge. A mirror or another person may assist with sizing.
2. The cushion opening will be directly under your nostrils. The cushion will hug your nose and mouth, and leak will be minimal.
3. At no time should your nose be inside of the cushion opening. If any part of the cushion overlaps your nose, the cushion is being worn incorrectly. A different size cushion may be needed.

Contact your healthcare provider for additional information.

---


### Philips Respironics DreamStation Mask Type and System One Resistance Control Settings

Using your mask with a Philips Respironics DreamStation or System One device provides optimal comfort. The provider sets this value (X1) on your device.

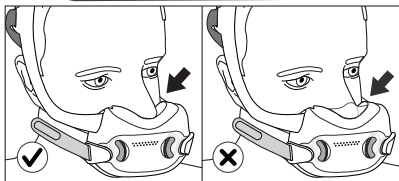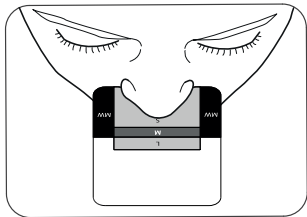
## Putting on the Mask

1. With the mask assembled (see **Assembly** section), place the cushion under the nose. ②
2. Position the frame on the top of your head. ③
- ✳️ **Tip:** For proper sizing, see **Cushion Sizing** and **Mask Frame Selection** sections.
3. Pull the headgear over the back of the head. ④
4. Attach the magnetic headgear clips to the mask cushion. ⑤
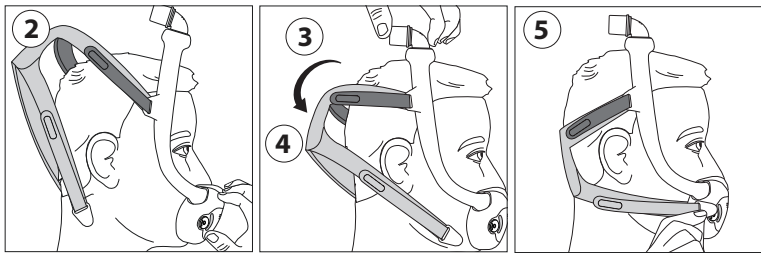
## Adjusting the Mask



1. Peel the headgear tabs away from the fabric. Adjust the length of the top straps. Press the tabs back against the fabric to reattach. Repeat this step with the bottom straps. ⑥
- ⚠️ **Note:** Do not overtighten the headgear. Signs of overtightening include redness, sores, or bulging skin around the edges of the mask.
2. Position the mask until it fits comfortably. ⑦
3. When finished the elbow should rest at the top of the head. ⑧

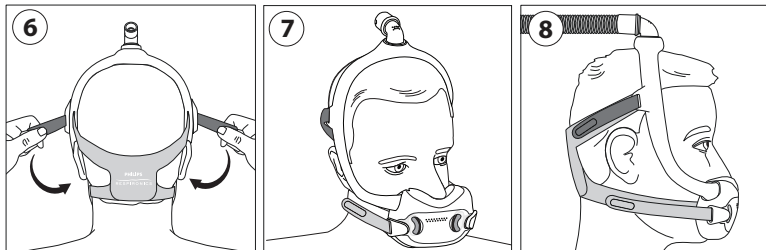
## Mask Frame Sizing

- Medium (MED) frame
- Small (SM) frame
- Large (LG) frame

The medium (MED) mask frame will comfortably fit most faces. If the MED frame does not fit your face, contact your healthcare provider to see if a small (SM) or large (LG) mask frame would better suit your needs.

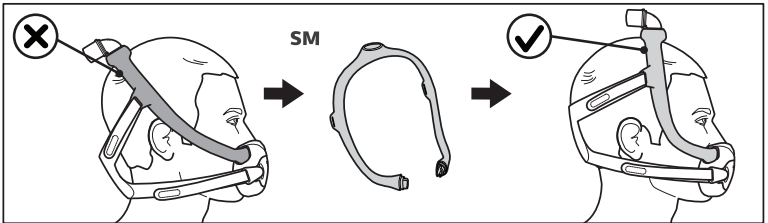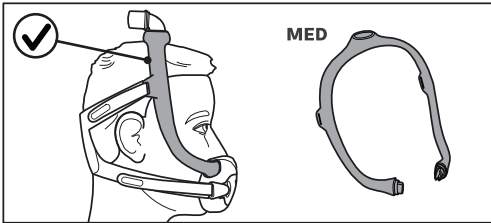✳️ **Tip:** If the mask frame falls backward and is too close to your ears you may need a smaller mask frame.

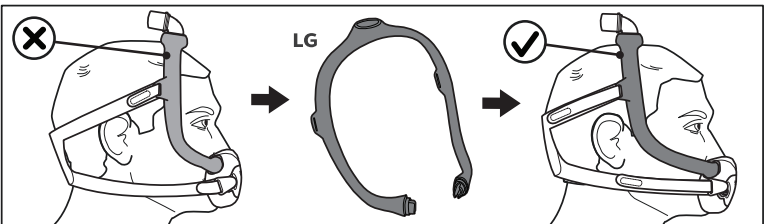✳ **Tip:** If the mask frame falls forward on the head and is too close to your eyes you may need a larger mask frame.

---

### Using the Mask

1. Connect the CPAP tubing (included with the CPAP device) to the tubing quick release on the elbow.
  2. Turn the therapy device on. Lie down. Breathe normally.
  3. Assume different sleeping positions. Move around until comfortable. If there are any excessive air leaks, make final adjustments. Some air leaking is normal.
- 


### Tube Quick Release

The elbow is equipped with a tubing quick release. Pull on the tubing quick release and it will disengage from the elbow and come loose with the tubing. ⑨

---

### Removing the Mask

To keep your adjustments, remove the mask by grasping the cushion and pulling forward away from the nose. Then pull the cushion and mask up and off of the head. ⑩
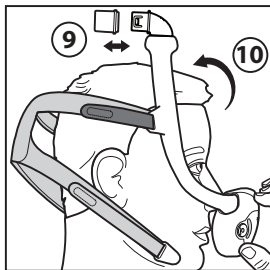
## 5 Mask Care

---

### Cleaning Instructions

**Daily:** Hand wash the non-fabric parts.

**Weekly:** Hand wash the fabric parts.

1. Disassemble the mask (see **Disassembly** section).
  2. Submerge and hand wash the mask in warm water with liquid dishwashing detergent.
- ⊖ **Note:** Ensure there are no air pockets present inside the mask parts while submerged.
3. Rinse thoroughly.
  4. Lay the headgear and fabric sleeves flat or line dry. Make sure the entire mask is dry before use.
- ⊖ **Note:** Do not place the fabric parts into a clothes dryer.
- ⚠ **Caution:** Any deviation from these instructions may impact the performance of the product.

**⚠ Caution:** Do not use bleach, alcohol, cleaning solutions containing bleach or alcohol, or cleaning solutions containing conditioners or moisturizers.

## Dishwashing Instructions

In addition to hand washing, the mask may be cleaned in the dishwasher once a week.

**⚠ Caution:** Use a mild liquid dish washing detergent only to wash the mask.

1. Remove the fabric parts. Do not wash the fabric parts in the dishwasher.
2. Wash in the top shelf of the dishwasher.

**⊖ Note:** Do not use the drying cycle on the dishwasher.

3. Air dry. Make sure the mask is dry before use.

**⚠ Caution:** Do not use bleach, alcohol, cleaning solutions containing bleach or alcohol, or cleaning solutions containing conditioners or moisturizers.

## 6 Assembly and Disassembly



### Assembly

1. Cushion: Press into the mask frame until the cushion clicks into place. **11****12**

**⊖ Note:** The end of the cushion and the opening in the mask frame are 'D' shaped. The cushion end should match the frame.

2. Elbow: Insert into the top of the mask frame.
3. Tubing quick release: Push the tubing quick release onto the elbow until it clicks into place.
4. (Optional) Fabric sleeves: Wrap the sleeves around the mask frame and press the strips together.

**⊖ Note:** The sleeve seams should be on the outside of the frame.

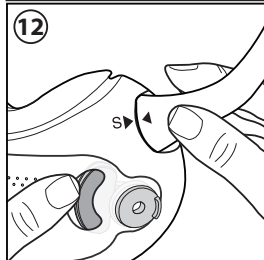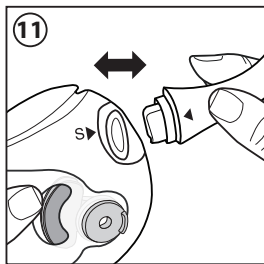
5. Headgear: Slide the headgear top strap tabs into the slots on the mask frame and fold backwards. Attach the magnetic headgear clips.

**Note:** The Philips Respironics logo will be on the outside and facing up when correctly assembled.

---

### Disassembly

1. Headgear: Undo the headgear tabs and pull through the slots in the mask frame. Turn the magnetic headgear clips away from the mask frame.
2. Fabric sleeves: Undo the strips and remove from the mask frame.
3. Tubing quick release: Pull the tubing quick release from the elbow.
4. Elbow: Pull the elbow from the top of the frame.
5. Cushion: Pull the cushion away from the mask frame.
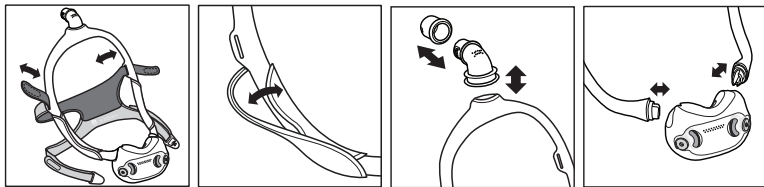
## 7 Healthcare Provider & Clinician Information

---

### Useful Life

The useful life of the DreamWear Full Face mask depends on the use conditions and maintenance (cleaning, institutional disinfection, and component replacement) of the device. Inspect the mask parts regularly for damage or wear. Replace components as deemed necessary.

---

### Multi-Patient Use

Clean and disinfect the mask daily and between patients when in use in a clinical setting. For multi-patient use in the hospital/institutional environment, use the Disinfection Guide to reprocess the mask between patients. These instructions can be obtained online at [www.philips.com/IFU](http://www.philips.com/IFU), or by contacting Philips Respironics Customer Service at 1-800-345-6443 (USA or Canada), at +1-724-387-4000, or go to [www.respironics.com](http://www.respironics.com) to find your local customer service contact information.

---


## Specifications

The technical specifications of the mask are provided for your healthcare professional to determine if it is compatible with your CPAP or bi-level therapy device.

**⚠ Warning:** Use outside of these specifications may result in ineffective therapy. The pressure flow curve shown below is an approximation of expected performance. Exact measurements may vary.

---

### Pressure Flow Curve
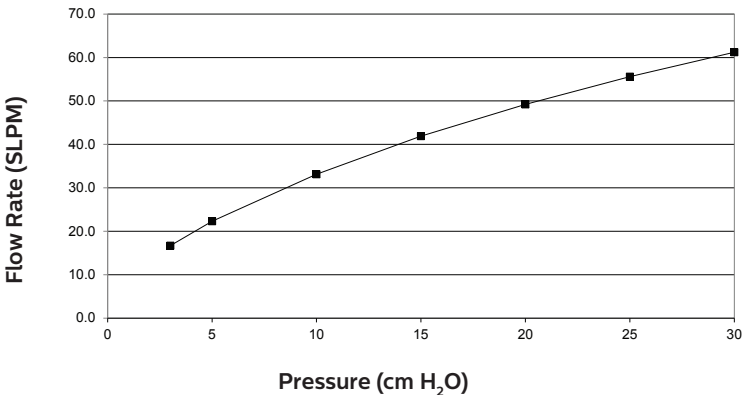
### Resistance with Flaps Closed to Atmosphere

Drop in Pressure (cm H <sub>2</sub> O)	at 50 SLPM	at 100 SLPM
All Cushion Sizes	0.9	3.1

---

**Closed to Atmosphere:** A minimum of 3 cm H<sub>2</sub>O must be maintained to ensure that the flaps close.

**Open to Atmosphere:** At an airway pressure less than -1 cm H<sub>2</sub>O, the flaps will be open.

---

### Inspiratory and Expiratory Resistance with Flaps Open to Atmosphere at 50 SLPM (cm H<sub>2</sub>O)

Inspiratory Resistance      0.4

Expiratory Resistance      1.8

---

---

## Deadspace (mL)

Frame size:	SM	MED	LG
Cushion size:	S: 176.3	S: 179.3	S: 182.2
	M: 182.7	M: 185.7	M: 188.6
	L: 201.3	L: 204.3	L: 207.2
	MW: 185.2	MW: 188.2	MW: 191.1

## Sound Levels (dBA)

A-weighted Sound Power Level 24.0

A-weighted Sound Pressure Level at 1 m 16.5

---

## Disposal

Dispose of the mask and headgear annually.

Dispose of this device in accordance with local collections and recycling regulations. For more information, visit [www.philips.com/recycling](http://www.philips.com/recycling).

---


## Storage Conditions

Temperature: -4° F to 140° F (-20° C to 60° C)

Relative Humidity: 15% to 95%, non-condensing
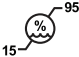
## 8 Symbols Glossary

Refer to <http://www.symbols.philips.com> for a description of the symbols used on this device and its packaging.

Symbol	Title and Meaning	Symbol	Title and Meaning
	Manufacturer - Indicates the medical device manufacturer		Date of manufacture - Indicates the date when the medical device was manufactured.  Country of manufacture - Indicates the country of manufacture of the product.  Note: When applied to the label, "CC" is replaced by the country code.
	Not Made with Natural Rubber Latex		
	Prescription device		
	Operator's manual; operating instructions Consult instructions for use.		Batch code - Indicates the manufacturer's batch code so that the batch or lot can be identified.
	Caution - Indicates the need for the user to consult the instructions for use.		Unique Device Identifier - Indicates the Unique Device Identifier information.
	Humidity limitation - Indicates the range of humidity to which the medical device can be safely exposed.		Temperature limit - Indicates the temperature limits to which the medical device can be safely exposed.
	Catalogue number - Indicates the manufacturer's catalogue number so that the medical device can be identified.		Philips Respironics DreamStation Mask Type and System One Resistance Control Value  Indicates the therapy device comfort setting.
	Packaging unit - Indicates the number of pieces in the package.		

## 9 Limited Warranty

---

Respironics, Inc., a Philips company (“Philips Respironics”) provides this non-transferable, limited warranty for DreamWear Full Face Mask (“Product”) to the customer who originally purchased the Product directly from Philips Respironics.

**What this Warranty Covers:** Philips Respironics warrants each new Product will be free from defects in materials and workmanship and will perform in accordance with the Product specifications under normal and proper use and maintenance in accordance with applicable instructions, subject to the exclusions below.

**How Long does this Warranty Last:** For the product, ninety (90) days from the longer of the date of shipment to the purchaser or date of setup by purchaser for the end user. The warranty period for accessories, replacement parts, and disposables including, but not limited to, chin straps is ninety (90) days from the date of shipment to the original purchaser.

**What this Warranty does not cover:** This warranty does not apply to any software included with the Product as the software warranty is included in the software license. This warranty does not cover damage or injury whether to the Products, personal property, or persons caused by accident, misuse, abuse, Acts of God, water ingress, repair or alteration by anyone other than Philips Respironics or its authorized service center, failure to operate in accordance with the terms of the operating manual and instructions, lack of reasonable care, the discontinuance of a network (e.g. 2G, 3G, etc.) by a carrier (e.g. ATT, Verizon, etc., or other defects not related to material or workmanship. This warranty is not transferable. If Philips Respironics finds that a Product returned for service or the issue raised is not covered under this limited warranty, Philips Respironics may charge an evaluation fee and return shipping.

**What Philips Respironics will do:** If a Product fails to conform to the warranties set forth above during the applicable warranty period, Philips Respironics will repair or replace the Product or refund the original purchase price, in Philips Respironics sole discretion. Philips Respironics may use new or remanufactured assemblies, components, and parts in repair and new or recertified refurbished devices for replacement. The balance of the original warranty period will apply to any Product or component of a Product repaired or replaced under this warranty.

**Warranty Disclaimer; Limitation of Liability:** EXCEPT AS SET FORTH IN THIS LIMITED WARRANTY, PHILIPS RESPIRONICS MAKES NO WARRANTIES, EXPRESSED OR IMPLIED, STATUTORY OR OTHERWISE, REGARDING THE PRODUCT OR ITS QUALITY OR PERFORMANCE. PHILIPS RESPIRONICS SPECIFICALLY DISCLAIMS THE IMPLIED WARRANTY OF MERCHANTABILITY AND THE IMPLIED WARRANTY OF FITNESS FOR A PARTICULAR PURPOSE. IN NO EVENT WILL PHILIPS RESPIRONICS MAXIMUM LIABILITY UNDER THESE WARRANTIES EXCEED THE ORIGINAL PURCHASE PRICE OR WILL PHILIPS RESPIRONICS BE LIABLE FOR ANY ECONOMIC LOSS, LOSS OF PROFITS, OVERHEAD, OR SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES.

Repair, replacement, or return of purchase price by Philips Respironics is the original purchaser's sole and exclusive remedy under this warranty.

This warranty gives you specific legal rights, and you may also have other rights that vary from state to state. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above exclusion and limitations may not apply to you.

**How to get warranty support:** Patients contact your local authorized Philips Respironics dealer and dealers contact Respironics, Inc. at:

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